


# Aortic Punch




Free-float™ cutting action  
ensures precise openings for every anastomosis.

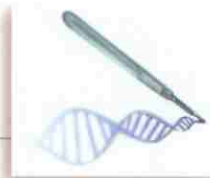

BRINGING  
BIOTECHNOLOGY  
TO SURGERY

**VENUS SURGICALS PVT. LTD.,**  
3362Y, AF - 8TH STREET, 11TH MAIN ROAD,  
ANNA NAGAR, CHENNAI - 600040.  
TEL : 91-44 - 26204748, 26203073, 26286956  
FAX : 91-44 - 26260665  
E-mail : [venussurgical@vsnl.com](mailto:venussurgical@vsnl.com)

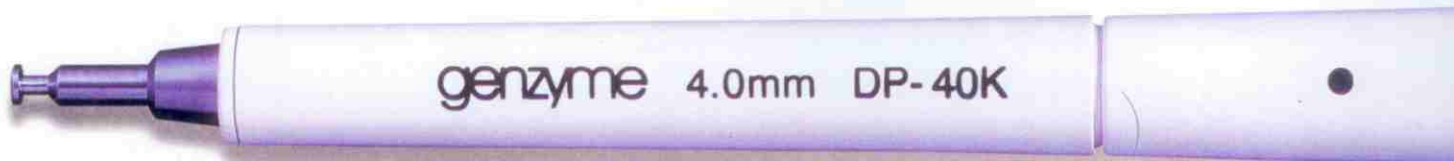

Clean, precise openings. Reduced frustration. Minimized costs. These can all be achieved by using the Genzyme Aortic Punch.

**Free-Float Cutting Action Technology**

By allowing the blade to "float" around the punch core during operation, alignment is maintained and mechanical jamming is prevented. One punch can be used to create multiple openings within a single procedure.

**Dual-Cutting Action**

Two diamond-sharp cutting surfaces lift and cleanly cut through all layers of tough tissue. The smooth cutting action of blades reduces trauma around the anastomotic site.

**Tapered Cutting Blade**

Tapering the distal portion of the cutting blade increases visibility around the anastomotic site.

**Increased Surgeon Comfort**

The Aortic Punch was designed with a generously sized thumb seat and finger grips. Distributing the cutting forces over a larger area of the hand greatly reduces discomfort.

**Functional Sizes**

Eight color-coded blade sizes create precisely the opening your patient needs. The standard length punch with a 7-3/4" handle enhances visibility. The mid-length punch with a 6-3/4" handle offers maximum maneuverability.

**Ordering Information**

For more information or to order the Aortic Punch that is right for your needs, contact your Genzyme Biosurgery representative or call 800-367-7874.

**Standard Length**

**Mid-Length**

**Standard Length Sterile**

Code Number	Size	Color	Quantity
DP-28K	2.8mm		Box of 6
DP-36K	3.6mm		Box of 6
DP-40K	4.0mm		Box of 6
DP-44K	4.4mm		Box of 6
DP-48K	4.8mm		Box of 6
DP-52K	5.2mm		Box of 6
DP-56K	5.6mm		Box of 6
DP-60K	6.0mm		Box of 6

**Mid-Length Sterile**

Code Number	Size	Color	Quantity
MDP-28K	2.8mm		Box of 6
MDP-36K	3.6mm		Box of 6
MDP-40K	4.0mm		Box of 6
MDP-44K	4.4mm		Box of 6
MDP-48K	4.8mm		Box of 6
MDP-52K	5.2mm		Box of 6
MDP-56K	5.6mm		Box of 6
MDP-60K	6.0mm		Box of 6

©2002 Genzyme Corporation  
All rights reserved. Printed in USA.  
Code No. 25-0863. 2002-04

**VENUS SURGICALS PVT. LTD.,**  
3362Y, AF - 8TH STREET, 11TH MAIN ROAD,  
ANNA NAGAR, CHENNAI - 600040.  
TEL : 91-44 - 26204748, 26203073, 26286956  
FAX : 91-44 - 26260665  
E-mail : venussurgical@vsnl.com



a division of  
GENZYME CORPORATION  
One Kendall Square  
Cambridge, MA 02139 USA  
800-367-7874  
617-252-7500  
www.genzymebiosurgery.com